

Selectable various attaching units for the therapeutic purpose.

■ Standard probe unit

● B probe unit



Irradiate to a deeper point

Output power : 1.9W
Focus area : 10mm

● C probe unit



Irradiate to a wide area

Output power : 2.5W
Focus area : 80mm

■ Option probe unit

● SG probe unit



Irradiate to a narrow area
such as the neck region

Output power : 1.3W
Focus area : 7mm

■ Product Specification

Category	: Infrared therapy equipment
Model name	: Superlizer mini PRO
Dimension	: W 290mm, D 250mm, H 120mm
Weight	: 1.8Kg
Electrical rating	: AC100V, 70VA, 50/60Hz, DC12V, 2.08A
AC adaptor	: AC100-240V 50/60Hz (Output:11-13Vdc 2.27-1.92A)
Operating temperature and humidity	: +10~+30℃, 30~70%RH
Rate time	: 10 minutes
Output power	: 2.5W



Before use for safety

- Read instruction manual before use.
- Please do maintenance of the equipment regularly.

Specification may be changed without notification.

■ Contact Information

 **Tokyoiken Co.,Ltd.** <https://www.tokyoiken.co.jp>

● Sales Department

2-27-2, Yushima, Bunkyo-ku, Tokyo 113-0034 Japan TEL +81-042-378-6630
9:00AM ~ 5:00PM Monday ~ Friday (Except Japanese holiday)

● Headquarter

1131-1, Higashi-Nagnuma, Ingaishi, Tokyo 206-0802 Japan TEL +81-042-378-6630



Compact Size & Lightweight Infrared therapy equipment




Excellent therapeutic effect!

—SUPER LIZER—
mini PRO

Superlizer mini PRO is a compact and lightweight infrared therapy equipment and specialized for the manual operated hand held therapy.

■ Beam irradiation therapy with Superlizer.

"Superlizer" is the phototherapy device that can irradiate a high power near-infrared light (0.6-1.6 μ m), a deep penetrating wavelength range, in spot. Its effects are recognized by various departments including pain clinics, and the device has attracted large attention from the medical community. The device meets various needs in the medical community such as operability, therapeutic effects, safety at high level.

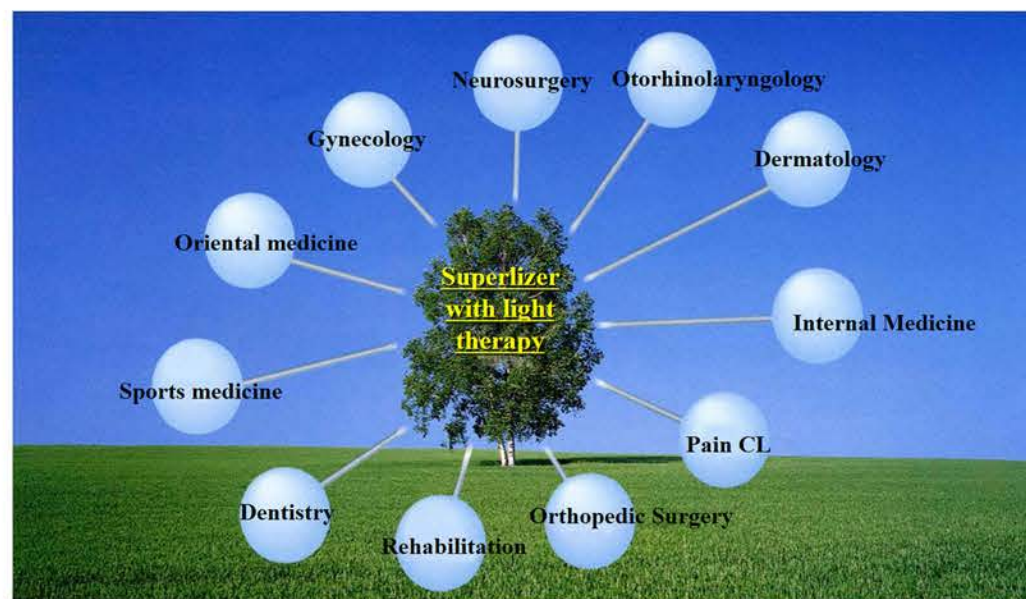


- Pain of Inflammation, An External Wound (Chronic Stage) Sprain, Fracture, Tendonitis, etc...
- Arthritis, Rheumatism Rheumatoid Arthritis, Shoulder Periarthritis, etc...
- Skin Disease Chronic Dermatitis, Atopic Dermatitis, Eczema, Pressure Ulcer, Skin Ulcer, Herpes Zoster, etc...
- Oral Surgery Temporomandibular Disorders, Masticatory Muscle Syndrome, etc...
- Others Raynaud's Disease, Whiplash, Myotonic Headache, Stiff Shoulder, Knee Arthropathy, Alopecia Areata, Sudden Deafness, Facial Paralysis, Autonomic Nerve Disease, etc...

■ Example of irradiation



1 Deep penetration into the human body by near infrared by 2,500mW maximum output power.

Extraction of near infrared (0.6 μ m~1.6 μ m) is deep penetration into the human body by the specialized optical filter, 2,500mW maximum output power into human body deeply makes warm, comfortable feelings and excellent therapy effect.

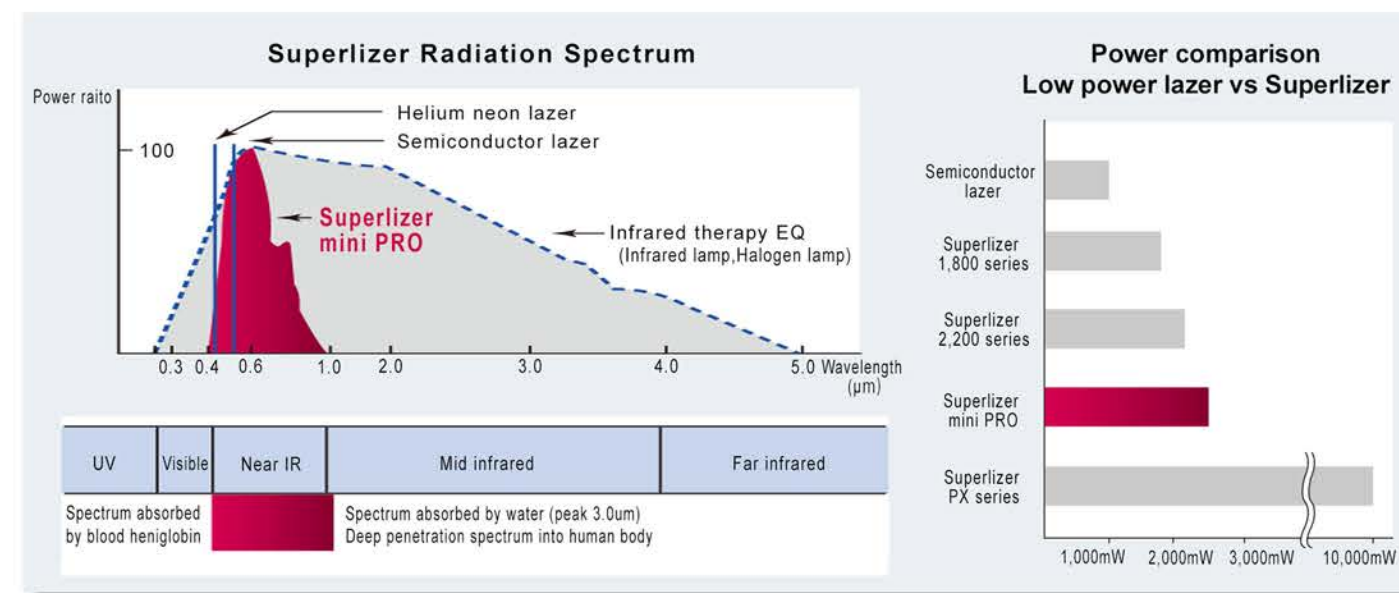

2 Simple operating panel. Two selectable irradiation modes.

Excellent therapy by combination of continuous mode & cycle irradiation mode. The device can be easily operated with dials.



3 High irradiation power by compact hand probe with internal halogen lamp.

Irradiation to various therapy points with 3 types of attaching probe units. Irradiate to a deeper point unit, wide area, narrow area such as the neck region.